



# ALOVA™

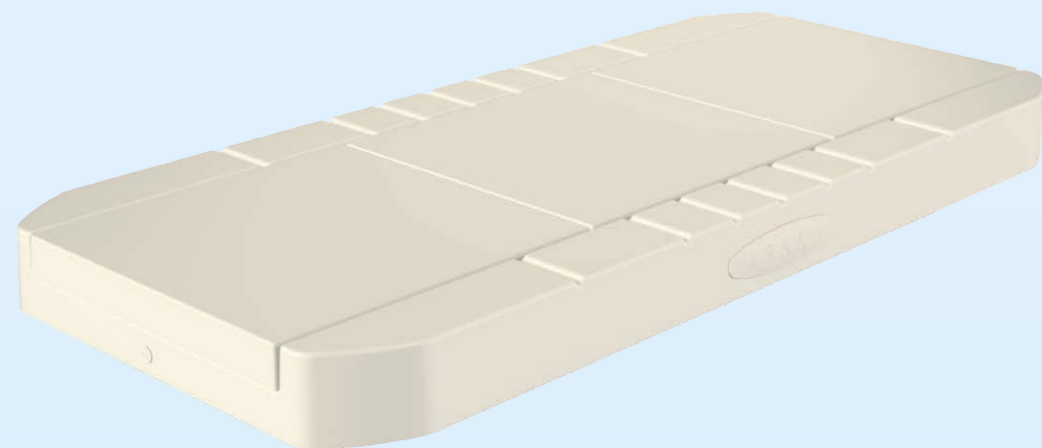

### Braden scale score

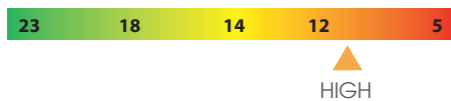

## KEY FEATURES

- ▶ Single-piece mattress with different hardnesses and stabilizing lateral edges made up of:
  - Molded visco elastic foam top layer: density:  $\geq 75 \text{ Kg/m}^3$
  - High Resilient (HR) polyurethane foam base: density:  $\geq 34 \text{ Kg/m}^3$ ; hardness: 3.5 kPa
- ▶ 3 different hardnesses: medium (1.7 kPa), firm (1.9 kPa) and soft (1.5 kPa)
- ▶ Patient weight: 30 to 135 Kg
- ▶ Dimensions: 195 x 88 x 14 cm (76.7 x 34.6 x 5.5 inches)
- ▶ Mattress weight: 19 Kg
- ▶ Covers available: Promust PU and PU HD

▶ Different hardnesses



## KEY BENEFITS

- ▶ Visco elastic foam provides excellent pressure re-distribution properties: homogeneous immersion of the body into the foam, thus increasing the bearing surface and reducing pressure applied in areas at risk of pressure ulcers.
- ▶ Different hardnesses adapted to each body area.
- ▶ The stabilizing lateral edges helps getting out of bed and patient transfer.
- ▶ Visco elastic foam provides patient comfort.

## INDICATIONS

HELP IN PREVENTING PRESSURE ULCERS FOR PATIENT AT MEDIUM TO HIGH RISK OF PRESSURE ULCER DEVELOPMENT, BEDRIDDEN MORE THAN 15 HOURS A DAY.

HELP IN TREATING EXISTING PRESSURE ULCERS FROM STADE 1 TO 2 COMBINED WITH TECHNICAL AID DEVICES.

PATIENT WEIGHT: 30 TO 135 KG/ 4.7 TO 21.2 ST

- ▶ Single-piece mattress made up of molded visco elastic foam top layer on High Resilient (HR) polyurethane foam base
- ▶ Stabilizing lateral edges
- ▶ Removable and washable integral cover



DESCRIPTION	ITEM REFERENCES ACCORDING TO PACKAGING CARDBOX	FIRE STANDARDS
ALOVA™ Promust PU cover	VMA/BS/HIP	EN 597-1 and 2 GPEM D1 90 Classe D
ALOVA™ Promust PU HD cover	VMA/BS/HIPH	

3 years warranty

Cleaning instructions: Promust PU HD cover



Photos not contractually binding

FP000036/01.2018

